

SWLFSYNERGYR™

PRESCRIPTION INFORMATION

SWLF Synergy R™ Prescription					0.018"	0.022"
Maxillary	Orientation	Angle	Torque	Offset	Part No.	Part No.
Central	R	3°	17°	0°	M07000	M07100
	L	3°	17°	0°	M07001	M07101
Lateral	R	9°	8°	0°	M07002	M07102
	L	9°	8°	0°	M07003	M07103
Cuspid	R	8°	0°	0°	M07006	M07106
	L	8°	0°	0°	M07007	M07107
1st/2nd Bicuspid	R/L	0°	-7°	0°	M07010	M07110
1st/2nd Bicuspid w/gingival offset	R/L	0°	-7°	0°	M07052	M07152
Mandibular	Orientation	Angle	Torque	Offset	Part No.	Part No.
Anterior (C/L)	R/L	0°	0°	0°	M07020	M07120
Anterior (Option)	R/L	0°	-6°	0°	M07022	M07122
Cuspid	R/L	0°	0°	0°	M07026	M07126
1st/2nd Bicuspid	R/L	0°	-17°	0°	M07030	M07130
1st/2nd Bicuspid w/gingival offset	R/L	0°	-17°	0°	M07062	M07162
Kits: 5-5 U&L with hooks on cuspids and bicuspids					KM07002	KM07102



1.800.525.6375

*The World's Oldest Synergistic, Bioprogessive,[®]
Breathing Enhancement Orthodontic Company.™*

SWLFSYNERGYR™



by rocky mountain orthodontics

SWLF SYNERGY R™

THE BEST JUST GOT BETTER


RMO's SWLF (Straight Wire Low Friction) Synergy R™ bracket represents the latest development in ConverTechnology: passive when you want it, total control when you need it. No clips, no doors, and no failures. SWLF Synergy R combines the simplicity and ease of self-ligating bracket design with the flexibility and advanced performance of Synergy's Friction Selection Control (FSC) modes. SWLF Synergy R provides minimal friction and rapid wire changeout, with cuspid and bicuspid brackets that can be converted into traditional Synergy® style brackets at any time during treatment. Clinically tested and proven effective, SWLF Synergy R is designed, engineered, and manufactured with pride in the USA.

Features and benefits include:

- cuspid and bicuspid brackets feature an integrated convertible cap
- can reduce treatment time and appointment intervals
- no moving parts—no broken clips, doors, or slides
- large flared lead-ins reduce kinking and binding
- low profile—comfortable for your patient
- convert to a standard Synergy style bracket at any time for advanced FSC modes
- consistent cap conversion strength—no undesired conversions
- rounded underwing clearance promotes easy ligation when converted
- permanent ID markings with vertical placement line for accurate positioning
- Metal-Injection-Molding (MIM) for consistent tolerances and smooth corners



Converting SWLF Synergy R bracket to Synergy style bracket with cap conversion pliers.

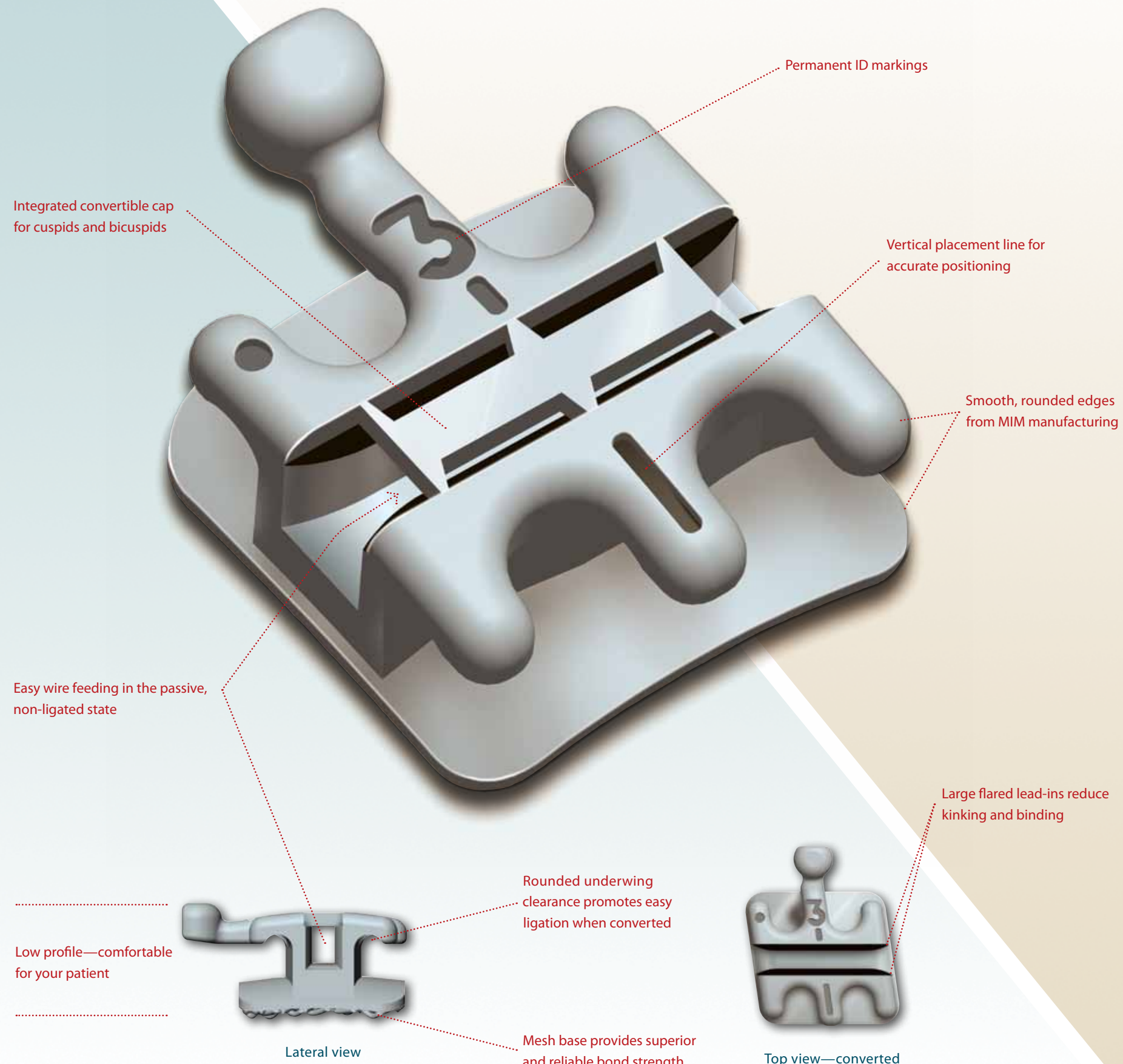

FRICITION SELECTION CONTROL

Take control of your treatment with FSC. Combined with SWLF Synergy R™'s integrated convertible cap, FSC modes deliver maximum tooth-by-tooth control throughout the entire course of treatment. Plus, clinicians can still satisfy color requests even during unconverted bracket stages by ligating the center wings without compromising performance. (Ligatures illustrated using original Synergy® bracket.)



REDUCED FRICTION



MODERATE ROTATION



MAXIMUM ROTATION



CONVENTIONAL CONTROL



MAXIMUM CONTROL